

# ANIOSRUB 85 NPC



Disinfecting solution  
for hydroalcoholic  
hand-rubbing



## • Concept SES :

### SAFETY

- ANIOSRUB 85 NPC is formulated without fragrance nor colouring.

### EFFICIENCY

- Wide antimicrobial spectrum: active against bacteria, yeasts, moulds and viruses.
- Hydroalcoholic solution to be used without water.

### SIMPLICITY

- Supports made of PVC specially made for hydroalcoholic products.

## INDICATIONS

Hydroalcoholic solution for hygienic treatment and surgical disinfection by hand-rubbing.

### Carry out hygienic treatment in these cases:

Removal of gloves. At the start of service/ at the end of service. A day-to-day activity, hotel activities. Treatments involving contact with healthy skin. After contact with a patient in septic isolation. Before carrying out any invasive procedure (peripheral catheter, urinary probe and other similar devices). After contact with infected patient or his/her environment. Between patients after any potentially contaminating procedure. Before any contact with a patient in protective isolation. Before carrying out a lumbar puncture, ascites puncture or articular drain-off, or other similar procedures. Before handling intravascular devices, pleural drains, implantable chambers and other similar situations. In the case of a series of contaminating actions for the same patient.

### Carry out surgical disinfection in these cases:

Before any surgical, obstetric or interventional radiology procedure. Before any procedure for which surgical asepsis is required: inserting a central or spinal catheter, implantable chamber, amniocentesis, pleural drain and other similar situations.

READY TO USE



## CHARACTERISTICS


- Ethyl alcohol: 755 ml/l
- Density at +20°C: approx. 0.870
- pH at +20°C: approx. 5.5 (near to skin pH)
- Aerobiocontamination study available on request.

# ANIOSRUB 85 NPC


Disinfecting solution for hydroalcoholic hand-rubbing

## INSTRUCTIONS FOR USE

Rub according to the EN 1500 method. **Do not rinse. Do not dry.**

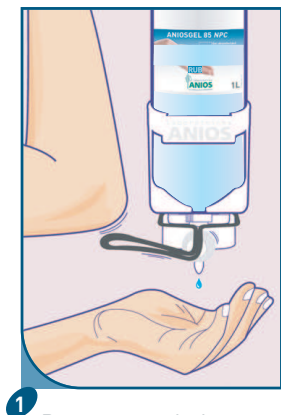

- 1 Pour as needed  
**Hygienic treatment:** sufficient amount to keep hands wet for 30 seconds  
**Surgical disinfection:** sufficient amount to keep hands wet for 2 x 45 seconds

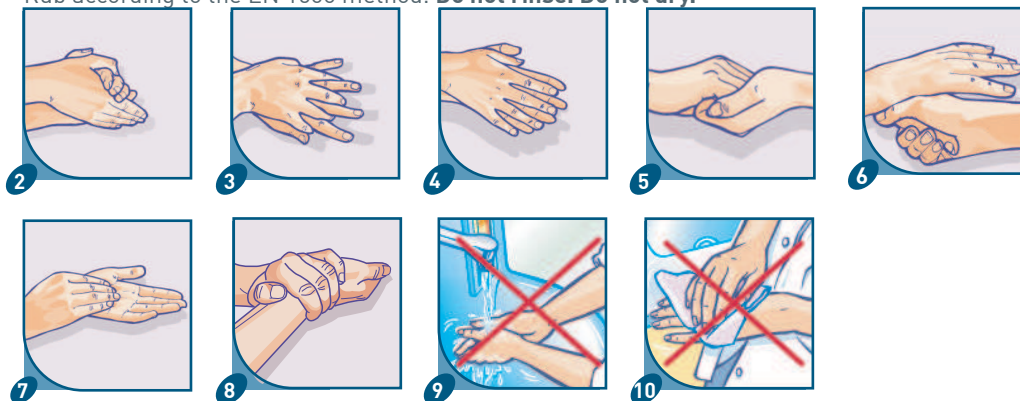

## COMPOSITION

Ethanol (700 mg/g ie 755 ml/l - N° CAS 64-17-5) in the presence of moisturizing and emollient agents, and water. Without perfume nor colouring.

## PRECAUTIONS FOR USE

Dangerous – respect the precautions for use (Drawn up according to the European rules in force regarding the classification and labelling of chemical products). Use biocidal products with caution. Before every use, read the label and information concerning the product. Store between +5°C and +25°C. Human hygiene biocidal product (Main Group 1 - PT1) – reserved for professional use.

## MICROBIOLOGICAL PROPERTIES

| Active against        | Standards                                                                                                                                                                                                                                                                                                        |
|-----------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Bacteria              | EN 1040, EN 13727 (2013) : 30 sec (hygienic friction) and 1 minute (surgical friction)                                                                                                                                                                                                                           |
| Mycobacteria          | EN 14348 : 30 sec                                                                                                                                                                                                                                                                                                |
| Yeasts / Moulds       | EN 1275, EN 13624 (2013) : 30 sec (hygienic friction) and 1 minute (surgical friction)                                                                                                                                                                                                                           |
| Viruses               | EN 14476 (2013) against Poliovirus : 2 min, Adenovirus and MNV : 30 sec<br>Activity on BVDV (surrogate of HCV), PRV (surrogate of HBV), HIV-1, HSV-1, Rotavirus, Vaccinia virus, RSV, Coronavirus (surrogate of SARS), Influenza virus (H <sub>1</sub> N <sub>1</sub> , H <sub>3</sub> N <sub>1</sub> ) : 30 sec |
| Hygienic disinfection | EN 1500 (2013) : Sufficient amount to keep hands damp for 30 seconds<br>• Report A13142-2 of 22/07/2013 – Hyggen Laboratory North - Prof. Med. A. Kramer<br>• Report SN 16982.1 of 19/12/2014 – Hyggen Laboratory - Prof. Med. H.-P. Werner                                                                      |
| Surgical disinfection | EN 12791 (2005): Sufficient amount to keep hands damp for 2 x 45 seconds<br>• Report A10060-2 of 11/03/2010 – Hyggen Laboratory North - Prof. Med. A. Kramer                                                                                                                                                     |

All the product's antimicrobial activity is included in the scientific dossier, available on request.

## PACKAGING

- 1 20 bottles, 75 ml each.....Ref. 1837.274
- 2 6 bottles, 300 ml each with screwed pump.....Ref. 1837.366
- 3 12 bottles, 500 ml each with screwed pump.....Ref. 1837.636\*
- 4 12 bottles, 500 ml each for dispenser.....Ref. 1837.639\*  
12 bottles, 1 L each for dispenser.....Ref. 1837.096
- 5 12 bottles, 1L each with screwed pump.....Ref. 1837.637\*
- 6 12 airless bottles, 1 L each.....Ref. 1837.638\*  
Supplier for 300 ml bottle.....Ref. 425.078  
Supplier for 500 ml bottle.....Ref. 425.058

\* Blue-coloured bottle



Pavé du Moulin  
59260 Lille-Hellemmes - France  
Tél. +33 3 20 67 67 67 - Fax : +33 3 20 67 67 68  
www.anios.com



GB1837-160212-Non contractual photographs